

S H O W
Y O U R
B O D Y
S O M E

Love

NEGATIVE PRESSURE
- POSITIVE RESULTS

 LYMPHATOUCH®

COMFORTABLE & NON-INVASIVE

LymphaTouch treatment reduces pain, swelling and bruising efficiently before and after surgery. It supports the lymphatic flow and blood circulation enhancing the natural healing process of the body. The unique combination of decompression and high frequency vibration is efficient also in treatment of fibrotic tissue and scars.

PAINLESS & RELAXING


The treatment is performed directly to the skin with various sizes of treatment cups to gently treat your patient from head to toe instead of an inconvenient treatment suit.

” LymphaTouch® has dramatically reduced tissue congestion, reduced swelling and inflammation, improved quality of life and reduced number of cellulitis.

Justine Whitaker

Director and Nurse Consultant of Northern Lymphology Ltd,
Senior Lecturer at the University of Central Lancashire, UK

More info available at lymphatouch.com

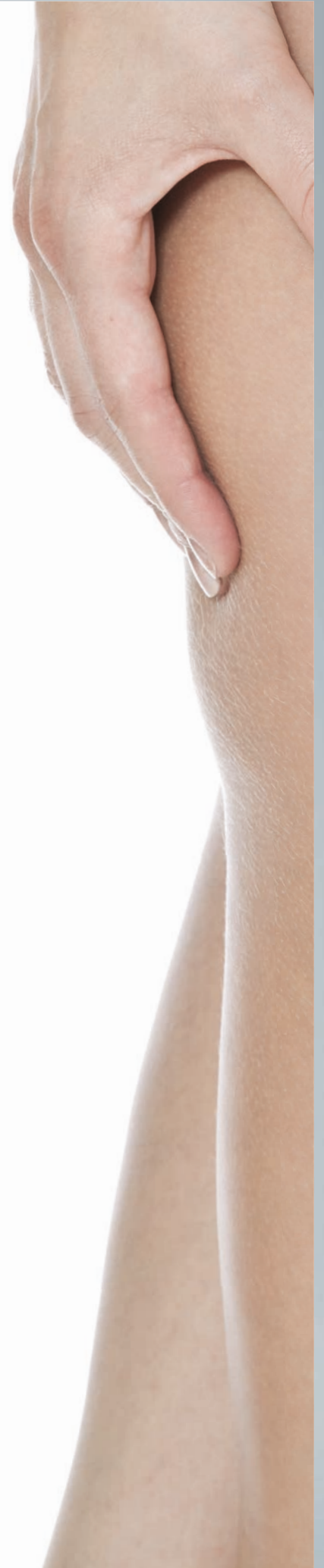



OPERATING PRINCIPLE

LymphaTouch® treatment makes a big difference for both professional and customer. Treatment areas include e.g. pre- and postoperative swelling and scarring, fibrosis, fascial tightness, active recovery, pain management and joint functionality.

LymphaTouch® expands and stretches the skin and the tissue pulling anchor filaments to dilate the endothelial openings of the lymphatic capillaries. Simultaneously, the expansion of fascial and connective tissue structures creates space for blood circulation and lymph flow.



” Investing in LymphaTouch® was definitely worth it. One year after purchasing the device my company has tripled its revenue and the number of new patients has increased significantly!

Leena Vallin

Nurse, pain management specialist, lymphtherapist and entrepreneur

3X MORE EFFICIENT THAN
MANUAL TREATMENT*

4X SOFTER SKIN THAN
MANUAL TREATMENT*

6X MORE AFFORDABLE**

1kg COMPACT SIZE
WITHOUT COMPROMISES

8H OF UNDISRUPTED USE

& NO TREATMENT OUTFIT

” LymphaTouch® is the most versatile product on the market. It can treat swelling, scar adhesion, and tissue restrictions caused by surgeries. It’s as easy as changing the setting on the screen.

Kathy McGinty
PTA, CLT, USA

*Supported by
clinical studies

** Comparison to
other leading
brands

ENHANCING RECOVERY



EFFICIENT IN SCAR TREATMENTS



CLINICALLY PROVEN IN VARIOUS PATIENT POPULATIONS



SUPPORTS THE BODY'S NATURAL HEALING PROCESS



NON-INVASIVE MECHANICAL DECOMPRESSION



RELIEVES PAIN, BRUISING AND DISCOMFORT



SAFE AND COMFORTABLE TO USE



#NaturalHealing

#NonInvasive

#Painless

#Safe

#BodyShaping

#Versatile

#ShowYourBodySomeLove

Find us on social media!



@LYMPHATOUCH

UNIQUE FEATURES OF LYMPHATOUCH®



TREATTOUCH 100®

Ergonomic and versatile treatment handle



DURABAT

Power supply for all day's use



APC

Adjustable pulsation and pressure control



INTELLITOUCH

Touchscreen UI with built-in treatment programs



TECHNICAL SPECIFICATIONS

DIMENSIONS

19 x 10 x 11,5 cm

DISPLAY

TFT LCD: 4.3" 800 x 480

POWER REQUIREMENT

Adapter: 100-240 VAC, 50-60 Hz

Output: 15Vdc, 2A, 30 W

TREATMENT CUPS

10 mm, 35 mm, 50 mm,

60 mm, 80 mm

STANDARDS

ISO 13485:2016

IEC 60601-1

IEC 60601-1-11

BATTERY

Rechargeable LT-1200 Lithium Ion

10,8 V, 5.2 Ah, 56 Wh

WEIGHT

1.1 kg

HIGH-FREQUENCY VIBRATION

Vibration range 20-90 Hz

NEGATIVE PRESSURE SETTING

Pressure range 20-250 mmHg,

possibility to upgrade to 350 mmHg

TREATMENT CUP FILTERS

Single use medical grade

air filters for treatment cups

CLASSIFICATION

Medical device CE-MDD/93/

42/EEC, Class 2a, CSA C22.2,

FDA class I Exempt, NRTL cert.

WARRANTY

2 years, excluding

batteries and accessories

